

# identibac

## *S. aureus* Genotyping Kit

Array Hybridisation Kit for DNA-based detection of resistance genes and pathogenicity markers of *Staphylococcus aureus* and assignment of unknown *S. aureus* isolates to known strains

***For Research Use Only. Not for Use in Diagnostic Procedures.***

## **CONTENT**

BACKGROUND .....	1
GENERAL INSTRUCTIONS FOR USE .....	2
Intended Use .....	2
Specifications .....	2
Technical Support .....	2
Safety Precautions .....	3
Material Safety Data Sheets (MSDS) .....	3
Shipping Precautions .....	3
REAGENTS AND DEVICES .....	4
Kit Components, Storage and Stability .....	4
Cell Lysis .....	4
DNA Labelling and Amplification .....	4
Hybridisation and Detection .....	5
Instrumentation & Software .....	5
Components Required but not Provided .....	6
PROTOCOL .....	8
Culturing and Harvesting Bacterial Cells .....	8
Extraction of DNA .....	8
Extraction of DNA by Spin Columns .....	9
Extraction of DNA by Automated Device .....	11
Linear Amplification and Internal Biotin Labelling .....	11
Hybridisation .....	12
General Remarks - Handling of Arrays .....	12
General Remarks - Handling of Liquids .....	13
General Remarks – the Substrate (Precipitating Dye) D1 .....	14
General Remarks - Thermoshakers .....	14
Protocol for Quantifoil's BioShake iQ .....	14
Adapted Protocol for Eppendorf's Thermomixer Comfort .....	16
Data Analysis .....	17
Starting the ArrayMate Reader .....	17
Worklist .....	18
Data Acquisition in the ArrayMate Reader .....	20
Results .....	22
Export of <i>S. aureus</i> Genotyping Test Reports .....	23

TROUBLESHOOTING .....	26
Staining Control.....	26
Image Quality .....	26
DNA Quality.....	27
Physical Damage to the Array .....	27
Ambiguous Results.....	27
Report Unavailable.....	28
Error Messages in Result Sheets.....	28
ADDITIONAL INFORMATION .....	30
Warranty .....	30
Disclaimer.....	30
Quality Control .....	31
List of Components for Separate Order.....	31
Legal Manufacturer.....	31
Contact .....	31
LITERATURE .....	32
UPDATES & SOFTWARE .....	32
APPENDIX 1 - FLOW CHART.....	33
APPENDIX 2 – IMAGES FOR TROUBLESHOOTING .....	35
APPENDIX 3 – TARGET GENES .....	36
APPENDIX 4 – TYPING INFORMATION .....	44
Definitions & Explanations.....	44
List of Currently Recognised Strains .....	45

## **BACKGROUND**

The CLONDIAG *S. aureus* Genotyping allows DNA-based detection of resistance genes and pathogenicity markers of *Staphylococcus aureus* and assignment of unknown *S. aureus* isolates to known strains.

RNA-free, unfragmented genomic DNA from pure and monoclonal *S. aureus* colony material is amplified approximately 40-fold and internally labelled with biotin-dUTP using a linear amplification protocol. In contrast to standard PCR, only one antisense primer per target is used resulting in producing single stranded (ss) DNA reaction products. This allows a simultaneous sequence specific labelling and amplification of an essentially unlimited number of targets. However, sensitivity is lower than in a standard PCR (whereas contamination with amplicons is nearly impossible) and for that reason the method is restricted to colony material and cannot be performed on samples such as swabs or pus. Resulting biotin labelled ssDNA is transferred and hybridised to DNA oligonucleotide microarrays with 333 probes for different genetic markers and a biotin staining control. These probes are printed in two spots each.

The target set consists of a variety of species markers, virulence-associated genes including genes that code for exotoxins, antibiotic resistance genes, genes encoding microbial surface components recognizing adhesive matrix molecules (MSCRAMMs), various enzymes and other types of markers [1].

Spot recognition is performed automatically based on a digital image of the arrays. The overall pattern is analysed automatically for the presence or absence of specific genes and it is compared to a database of strain profiles allowing assignment to clonal complexes and strains [1, 2].

## **GENERAL INSTRUCTIONS FOR USE**

### **Intended Use**

#### ***For Research Use Only. Not Intended for Use in Clinical Diagnostics.***

This kit allows genotypic characterisation of *S. aureus* isolates for research and epidemiological applications. It must not be used as a substitute for phenotypic susceptibility tests and for the guidance of antibiotic therapy. It cannot be used for other bacteria than *S. aureus*.

### **Specifications**

Upon receipt, the assay components need to be stored at different temperatures as specified on the package insert. The assay is to be performed at an ambient temperature of 18°C to 28°C.

### **Technical Support**

If you require any further information on this product please contact:

email: [cct.home@clondiag.com](mailto:cct.home@clondiag.com)

phone: +49 (0) 3641 3111 0

### **Safety Precautions**

*The kit is intended for use by personnel that are trained in microbiological and molecular methods. Preparation of DNA from pure *S. aureus* colonies (clones) requires expertise in microbiology and the local regulations for handling of pathogenic microorganisms (biosafety level 2) are to be obeyed.*

*Isolated, cell-free *S. aureus* DNA may be processed without further biosafety precautions, although contamination with *S. aureus* or other bacteria needs to be ruled out.*

*Always wear protective clothes as required for laboratory work by your local regulations.*

### **Material Safety Data Sheets (MSDS)**

According to OSHA 29CFR1910.1200, Commonwealth of Australia [NOHSC: 1005, 1008(1999)] and the latest amendments to the European Union Directives 67/548/EC and 1999/45/EC, the enclosed reagents do not require a Material Safety Data Sheet (MSDS). They do not contain more than 1% of a component classified as hazardous and do not contain more than 0.1% of a component classified as carcinogenic. MSDS therefore are not provided. Nevertheless, the buffers may cause irritation if they come into contact with eyes or skin, and may cause harm if swallowed. The regular precautions associated with laboratory work should be obeyed (e.g., wear protective goggles, gloves and lab coat and avoid contact with the reagents). In case, any liquids are spilled, clean with disinfectant and/or laboratory detergent and water.

Alere assumes no liability for damage resulting from handling or contact with these products. If you have any questions please contact our Technical Support (see above).

### **Shipping Precautions**

RID/ADR: *Kein Gefahrgut* / No dangerous goods

IMDG: No dangerous goods

IATA: No dangerous goods

## **REAGENTS AND DEVICES**

### **Kit Components, Storage and Stability**

All reagents are provided in a certain surplus amount (see below). In case of need, all components may also be ordered separately; please refer to the order numbers at the end of this manual. For pricing please contact your local representative or our customer service, respectively.

The expiry date can be found on each bottle and on the outer package. All components have been tested for stability for short term shipment (<1 week) at ambient temperature (< 37 °C). The kit components with a rather limited stability are D1 and C3. The other components have proven to be stable even six months after the kit expiry date has passed.

### **Cell Lysis**

- A1: Lysis Buffer  
Store at 18-28°C (ambient temperature). Surplus: 50%.
- A2: Lysis Enhancer (lyophilised)  
Store at 18-28°C (ambient temperature). Centrifuge A2 tubes shortly prior to opening. Add 200 µL Buffer A1 to Lysis Enhancer before use. Mix well and store for less than 1 week at 2-8°C. Sufficient for 96 isolations.

### **DNA Labelling and Amplification**

- B1<sup>ST</sup>: Labelling Buffer/Master Mix  
Store at 2-8°C. Surplus: 25%.
- B2: Labelling Enzyme  
Store at 2-8°C. Surplus: 50%.

## Hybridisation and Detection

- ArrayStrips (12 x 8 samples),  
Protected against light and sealed under inert gas. Store at 18°C to 28°C. After opening to be used within two weeks. Close the unused wells with caps, protect them against humidity and dust and store them at a dark place. *Avoid any touching or scratching of the surface of the microarray at the bottom of the well. Do not store or handle unused wells at an air humidity of more than 60% since this may irreversibly corrode the spots.*
- CapStrips (24 strips)
- C1: Hybridisation Buffer  
Store at 18-28 °C, protect against sunlight. Surplus: 100%.
- C2: Washing Buffer 1  
Store at 18 °C - 28 °C, protect against direct sunlight. Surplus: 100%.
- C3: HRP Conjugate 100x  
Store at 2-8 °C, protect against direct sunlight. Surplus: 100%.
- C4: Conjugate Buffer  
Store at 18°C to 28°C, protect against direct sunlight. Surplus: 200%.
- C5: Washing Buffer 2  
Store at 18°C to 28°C, protect against direct sunlight. Surplus: 200%.
- D1: Horseradish Peroxidase Substrate  
Store at 2-8°C, protect against direct sunlight. Surplus: 50%.
- CM: Reference DNA from *S. aureus* strain N315  
(GenBank accession number BA000018), 0.1 µg/µL. Store at 2-8 °C. Sufficient for 5-6 tests.

## Instrumentation & Software

- ArrayMate Reader (to be ordered separately, for details see below)  
The *S. aureus* Genotyping kit may be used on the ArrayMate reader only. The older devices



ATR01/03 are not suitable for reading ArrayStrip based assays. In case of any questions please contact your local distributor and/or Alere Jena.

- Iconoclust software (provided with the reader)
- Test specific software plug-in that contains information such as spot names, marker names, location of the spots on the array. This plug-in is delivered with the reader or can be downloaded from our website. Test-specific plug-ins will occasionally be updated. Please check the NEWS section of our website <http://www.clondiag.com>. Support is available under [cct.home@clondiag.com](mailto:cct.home@clondiag.com) or +49 (0) 3641 3111 0.

### **Components Required but not Provided**

- Growth media for the cultivation of *S. aureus*. The test should be performed with colonies harvested from Columbia Blood Agar. Other media that contain blood may also suffice, but have not systematically been tested. Media that do not contain blood (Mueller-Hinton, or MRSA selective media) usually yield lower DNA concentrations and should not be used. Liquid media should also not be used because contaminations or mixed cultures cannot easily be ruled out.
- Equipment and consumables needed for the cultivation of *S. aureus* (incubator, inoculation loops, Petri dishes)
- Clumping factor/coagulase assays for confirmation of *S. aureus*
- DNA preparation kit: The assay has been tested with the DNeasy Blood & Tissue Kit from Qiagen (cat# 69504), QIAamp Minikit (cat# 51306) and a DNA preparation kit for Qiagen's EZ1 automated device (DNA Tissue Kit, cat# 953034).  
Please note: DNA isolation from *S. aureus* requires a pre-treatment with the Cell Lysis components A1/A2 (see below).
- Equipment needed for DNA isolation, e.g. pipettes, centrifuge, thermoshaker or automated device (see above)
- Photometer for measuring the concentration of DNA

- Equipment for DNA gel electrophoresis for quality control of DNA
- Thermocycler
- Thermoshaker

We strongly recommend the BioShake iQ by Quantifoil Instruments (<http://www.qinstruments.com/>) equipped with a customised heating block designed to fit ArrayStrips. Alternatively, you may use Eppendorf's Thermomixer Comfort, equipped a heating block for microtitre plates.

- Pipettes: suitable for 1µL-5µL volumes, 90µL, 100µL, 200µL, 1000µL
- Multichannel Pipettes for 100-200 µL
- Reagent tubes suitable for PCR
- Ultrapure (PCR grade) water


## **PROTOCOL**

### **Culturing and Harvesting Bacterial Cells**

***S. aureus is a potential pathogen. All procedures for cultivation of the bacterium and DNA preparation need to be performed by properly trained staff in a biosafety level 2 facility.***

Grow *S. aureus* on Colombia blood agar (overnight at 37°C or 48 hrs at room temperature). Obtain confirmation of the identification as *S. aureus* (by katalase + coagulase/clumping factor assays or by other means) and make sure that you have a pure, monoclonal culture of *S. aureus*. Contamination with other bacteria, especially with other staphylococci needs to be strictly avoided as they might carry the same resistance genes as certain *S. aureus* strains and thus can introduce false positive signals and patterns.

- Centrifuge A2 tube shortly, open it, add 0.2 mL of Lysis Buffer A1 to Lysis Enhancer A2 and dissolve.
- Add an inoculating loop full of monoclonal colony material of the *S. aureus* isolate to this A1/A2 reagent, vortex.

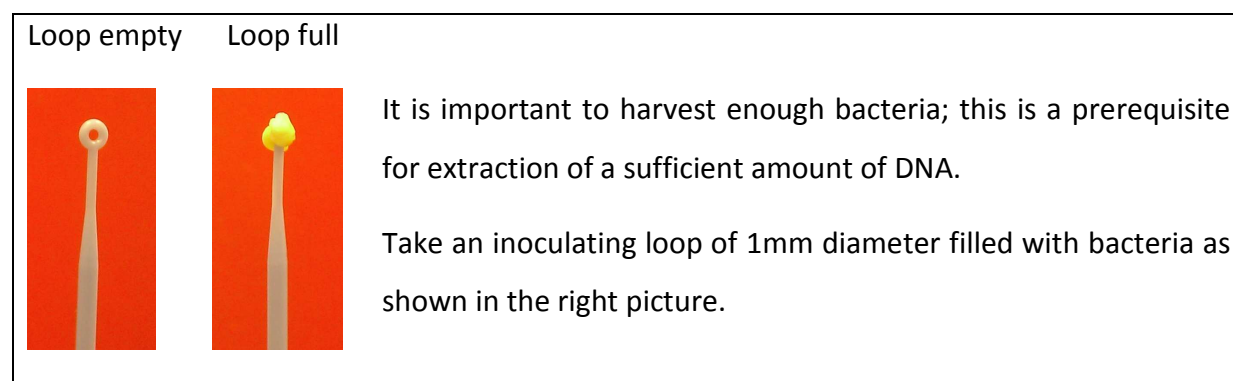

### **Extraction of DNA**

*The required sample type for the S. aureus Genotyping assay is 0.5-2 µg of intact genomic DNA from a single clone of S. aureus.*

This is much more DNA than for standard PCR applications (see Introduction).

*The DNA specimen needs to be free of RNA and it should not be fragmented.*

This can be determined by agarose gel electrophoresis. DNA should not be prepared by disrupting *S. aureus* cells using bead beaters, ultrasonication or aggressive chemicals such as in alkaline lysis protocols. Most performance problems with the *S. aureus* Genotyping kit are due to insufficient amounts or quality of DNA preparation. We therefore strongly recommend following the protocols outlined below.

We also recommend performing an experiment with Reference DNA from *S. aureus* strain N315 (CM reagent) when establishing the procedure in your lab. This will help to determine the cause of possible problems.

#### **Extraction of DNA by Spin Columns**

- Incubate the colony material of the *S. aureus* isolate in A1/A2 for 30-60 min at 37°C and 550 rpm in the thermoshaker.
- Proceed with the DNA preparation protocol of the DNA preparation kit. For the Qiagen DNeasy Blood&Tissue Kit that is as follows:
- Add 25 µL proteinase K (Qiagen Kit, or equivalent) and add 200 µL buffer AL (Qiagen Kit)
- Vortex shortly or shake vigorously.
- Incubate for 30-60 min at 56°C and 550 rpm in the thermoshaker.
- Add 200 µl ethanol (96-100%).
- Vortex the sample and centrifuge shortly.
- Transfer the complete content of the tube (including any precipitate) into a spin column that is placed in a 2 ml collection tube.
- Centrifuge at room temperature, time and speed need to be determined depending on viscosity of the sample and type of centrifuge used. All liquid should be collected in the collection tube afterwards.
- Discard collection tube with liquids.

- Place the spin column in a new 2 ml collection tube (provided with the kit).
- Add 500 µl Buffer AW1.
- Centrifuge at room temperature.
- Discard collection tube with liquids.
- Place the spin column in a new 2 ml collection tube (provided with the kit).
- Add 500 µl Buffer AW2.
- Centrifuge at room temperature, the membrane of the spin column should be dry, and all liquid should be in the collection tube.
- Discard collection tube with liquids.
- Place the spin column in a clean 1.5 ml tube (provided with the kit).
- Add 50 µl Buffer AE (or PCR grade distilled water) directly onto the membrane of the spin column.
- Incubate at room temperature for 5 min to elute DNA.
- Centrifuge.
- Optional: add another 50 µl Buffer AE (or PCR grade distilled water) directly onto the membrane, incubate at room temperature for 1 min and centrifuge again.
- Discard the spin column.

*Ethanol from Washing Buffers strongly inhibits the enzymes used in the assay.*

A contamination with Washing Buffer might occur during elution of prepared DNA by drops adhering to the funnel of the spin columns. Thus these funnels should be gently touched and tried with sterile filter paper or wipes prior to the elution step. Alternatively, prepared DNA can shortly be heated to evaporate ethanol (e.g., 10 min at 70°).

- Check for DNA integrity and absence of RNA (e.g., agarose gel). If necessary, you might perform another digestion step with additional RNase (not provided). Measure DNA concentration ( $A_{260}$  method), it shouldn't be less than 0.1 µg/µl. The concentration might be increased by heating and evaporation of water, or by using a speed vac centrifuge.

### **Extraction of DNA by Automated Device**

The assay has been tested with Qiagen's EZ1. Other systems also can be used. However, performance should be checked with some known reference strains prior to routine use. Incubate the colony material of the *S. aureus* isolate in A1/A2 for 30-60 min at 37°C and 550 rpm in the thermoshaker as described above (depending on the input sample volume required by the device you are actually using, the A1/A2 mixture might be divided into two aliquots, and used for DNA preparation of two samples).

- Add 10 µl proteinase K and add 100 µl buffer AL.
- Vortex shortly or shake vigorously.
- Incubate sample, 45-60 min at 56 °C and 550 rpm in the thermomixer.
- When the cells are lysed, proceed by performing the tissue lysis protocol (Bacteriacard) for Qiagen's EZ1
- *For Qiagen's EZ1:* Front row: empty elution tubes (1.5 ml); Second row: tip holder with tips; Third row: empty; Back row: sample tube with conical tip (2 ml) with the 200 µl sample volume. Set tissue lysis protocol with a set sample volume of 200 µl and an elution volume of 50 µl.
- Concentrate DNA and evaporate traces of solvents by heating the sample, at 70 °C for 5-10 minutes.

### **Linear Amplification and Internal Biotin Labelling**

Please keep in mind the limited surplus of reagents whilst pipetting. The surplus of B1<sup>ST</sup> labelling reagent is 25%.

- Prepare a Master Mix by combining 4.9 µL of B1<sup>ST</sup> labelling reagent and 0.1 µL of B2 (DNA polymerase) per sample.
- Add 5 µl DNA (0.5-2 µg) prepared as described above to a 5 µL aliquot of the mastermix. Do not forget to label the vial!

- Perform amplification in a pre-programmed thermocycler (such as Mastercycler gradient with heated lid, VWR, cat No. 460-0108) according to following protocol:

Pre-heat cover/lid to 105°C	
300 sec at 96°C	
50 (to 55) cycles with:	60 sec at 96°C
	20 sec at 50°C
	40 sec at 72°C
Cool down to 4°C, hold	

*Please note: When using another device, some adaptations, such as an increase of the number of cycles, might be necessary. Before starting routine use, please test the protocol with a few known reference strains and the control DNA supplied with the kit.*

## **Hybridisation**

### **General Remarks - Handling of Arrays**

*Never touch the array surface !*

*Avoid complete drying of the array surface during processing !*

*Do not allow it to stay without liquid for more than two minutes !*

*Never rinse the wells with distilled water after hybridisation.*

Unused wells should be capped during the whole procedure. The strips may be processed up to three times without a loss of quality of properly capped unused arrays. Close all wells that will not be used with a cap und leave it there until you use these wells (for storage conditions after use: see section “Kit components, storage and stability/Hybridisation and Detection”).

Always label your array strips with a laboratory marker at the recommended position. Never label them on the bottom or across the data matrix barcode! This may cause an error.




Avoid contact of data matrix barcode with organic solvents! The ArrayMate needs the information encoded in the data matrix to perform the assay.

Avoid touching the bottom of the microarray strip and keep it clean.


### General Remarks - Handling of Liquids

We recommend the use of a multichannel pipette and reagent reservoirs. Please keep in mind the limited surplus of C1 (100%).

We strongly recommend that the liquid is removed by pipetting rather than by inverting the strips and flicking the liquids out. Fine tipped soft, disposable Pasteur pipettes are suited best (such as VWR Best.nr: 612-2856).

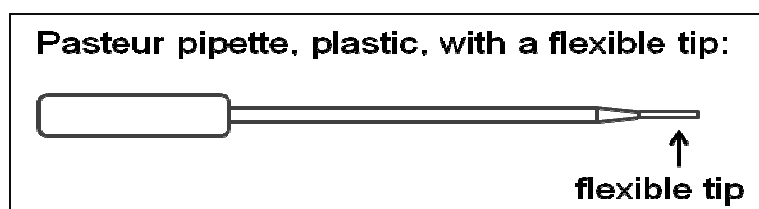

Always place the pipette tip at the cavity between the array and the wall of the reagent well. If you touch the array surface, probes may be scratched off and this may cause an error.

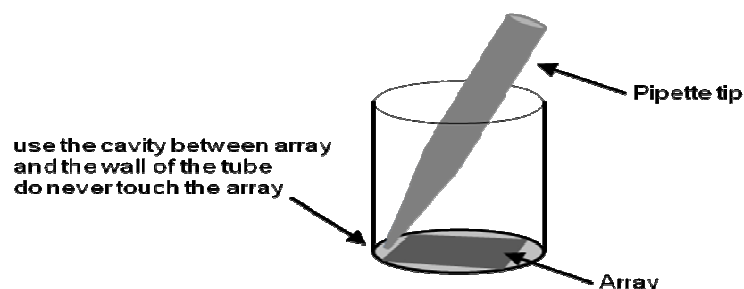



### **General Remarks – the Substrate (Precipitating Dye) D1**

An appropriate amount of substrate (precipitating dye) should be filled into an Eppendorf tube and taken out of the refrigerator when starting the procedure allowing it to pre-warm to room temperature/25°C. Cold D1 may yield weak signals. D1 should shortly be centrifuged prior to use to remove bubbles as well as possible precipitates.

Triggered by peroxidase, in case of positive reactions, the dye precipitates but it is not covalently bound. The precipitate can be dissolved by vigorous shaking. Thus the arrays must not be shaken, dropped or moved abruptly during the staining procedure and afterwards.

After completion of staining, remove and discard reagent D1 as completely as possible and scan immediately. The dye precipitate fades slowly in presence of liquids.

### **General Remarks - Thermoshakers**

The correct temperature within the vessels is essential; therefore always use the appropriate equipment for heating. Because of a possibly inhomogeneous distribution of the temperature within the heating block and because of possible differences between displayed and actual temperatures, the use of different brands of thermoshakers might affect test performance. We tested the assay with BioShake iQ by Quantifoil Instruments (<http://www.qinstruments.com/>) equipped with a customised heating block designed to fit ArrayStrips and Eppendorf's Thermomixer Comfort, equipped with a heating block for microtitre plates. Thus we recommend the use of either device. Accordingly, two slightly different protocols are discussed here.

When using other devices, some modifications to the protocol might be necessary. Before starting routine use, please test the protocol with a few known reference strains or the control DNA supplied with the kit.

### **Protocol for Quantifoil's BioShake iQ**

- Switch on the thermoshaker and let it pre-heat to 50 °C.
- Remove the ArrayStrip(s) from the pouch.

- Insert the ArrayStrip(s) into the white frame. Assure the correct orientation (data matrix barcode close to row (A) and proper fit.
- Pre-wash the array in two steps:
- First, PCR-grade distilled water, 150 µl per well at 50°C, 5 min and 550 rpm
- Second, C1 Hybridisation Buffer, 150 µl per well at 50°C, 5 min and 550 rpm
- Add 90 µL of buffer C1 to each tube with (10 µL) labelled amplification product, mix gently
- Remove the buffer from the well with the array and add the mixture of C1 and labelled amplification product
- Incubate at 50°C, 60 min and 550 rpm.
- Remove liquid and wash with 150 µl C2 Washing Buffer, pipett up and down four times, remove and discard.
- Add another 150 µl C2 Washing Buffer. Incubate at 30°C, 10min and 550rpm.
- Meanwhile, prepare conjugate: for each experiment add 1 µl conjugate 100xHRP to 100 µl C4 Conjugation Buffer. This mixture is stable for one day at room temperature; C3 is delivered with a surplus of 100%, C4 is delivered with a surplus of 200%.

Suggested pipetting scheme:

	1 well	2-3 wells	4-6 wells	7-10 wells	11-15 wells	16-20 wells	21-30 wells	31-40 wells
C3	1.5 µL	3.5 µL	7 µL	11 µL	16 µL	21 µL	32 µL	42 µL
C4	150 µL	350 µL	700 µL	1100 µL	1600 µL	2100 µL	3200 µL	4200 µL

- Remove and discard Washing Buffer, and add 100 µl diluted conjugate to each well, incubate at 30°C, 10min and 550rpm.
- Remove liquid and wash with 150 µl C5 Washing Buffer, pipett up and down four times, remove and discard.
- Add another 150 µl C5 Washing Buffer. Incubate at 30°C, 2 min and 550rpm.

- Remove and discard Washing Buffer, add 100 µl of D1 substrate (precipitating dye, at 25°C, see above) per well.
- Incubate at 25°C, 6 min *but do not shake* !
- Remove liquid completely.
- The bottom of the ArrayStrips may be cautiously be cleaned with wipes, bubbles may be removed by removing and adding D1.
- Scan and process (see below).

#### **Adapted Protocol for Eppendorf's Thermomixer Comfort**

- Switch on the thermoshaker and pre-heat to 55 °C.
- Remove the ArrayStrip(s) from the pouch.
- Insert the ArrayStrip(s) into the white frame. Assure the correct orientation (data matrix barcode close to row (A) and proper fit.
- Pre-wash the array in two steps:
- First, PCR-grade distilled water, 150 µl per well at 55°C, 5 min and 550 rpm
- Second, C1 Hybridisation Buffer, 150 µl per well at 55°C, 5 min and 550 rpm
- Add 90 µL of Buffer C1 to each tube with (10 µL) labelled amplification product, mix gently.
- Remove the buffer from the well with the array and add the mixture of C1 and labelled amplification product.
- Incubate at 55°C, 60 min and 550 rpm.
- Remove liquid and wash with 150 µl C2 Washing Buffer, pipett up and down four times, remove and discard.
- Add another 150 µl C2 Washing Buffer. Incubate at 30°C, 5min and 550rpm.
- Meanwhile, prepare conjugate: for each experiment add 1 µl conjugate 100xHRP to 100 µl C4 Conjugation Buffer. This mixture is stable for one day at room temperature; C3 is delivered with a surplus of 100%, C4 is delivered with a surplus of 200%.

Suggested pipetting scheme:

	1 well	2-3 wells	4-6 wells	7-10 wells	11-15 wells	16-20 wells	21-30 wells	31-40 wells
C3	1.5 µL	3.5 µL	7 µL	11 µL	16 µL	21 µL	32 µL	42 µL
C4	150 µL	350 µL	700 µL	1100 µL	1600 µL	2100 µL	3200 µL	4200 µL

- Remove and discard Washing Buffer, and add 100 µl diluted conjugate to each well, incubate at 30°C, 15min and 550rpm.
- Remove liquid and wash with 150 µl C5 Washing Buffer, pipette up and down four times, remove and discard.
- Add another 150 µl C5 Washing Buffer. Incubate at 30°C, 2 min and 550rpm.
- Remove and discard Washing Buffer, add 100 µl of D1 substrate (precipitating dye, at 25°C, see above) per well.
- Incubate at 25°C, 6 min *but do not shake* !
- Remove liquid completely.
- The bottom of the ArrayStrips may be cautiously be cleaned with wipes, bubbles may be removed by removing and adding D1.
- Scan and process (see below).

## **Data Analysis**

### **Starting the ArrayMate Reader**

We recommend starting the ArrayMate Reader after having started the hybridisation; this allows you to conveniently start the device and to import the worklist file (see below).

Please note that this is a short instruction only. For more detailed information please refer to the ArrayMate User Manual.

- Switch on the ArrayMate (main switch on the rear below the electric cable plug, operating switch on the bottom/left corner of the front side).
- Switch on the screen (switch right hand below the screen).
- Log-in as “R&D User” (Research and Development User) for full access to test specific software (a default password will be provided together with the ArrayMate device). If you log-in as “User”, you will obtain only raw values, but no interpretation as positives/negatives and no strain assignment. “Administrator” log-in will allow manipulation of file folders and software; and this should be done only upon direct advice of Clondiag’s IT team.
- The user interface will be loaded, ArrayMate performs internal testing. This will require slightly less than a minute.
- Click on the icon “New Run” (left upper edge of the screen). A suggestion for a run name / folder name for the new run appears in the top line of the screen). You may modify or change the experiment name at your convenience.
- Type in your operator ID (optional).
- You may enter a comment into the “memo” field (optional).

## Worklist

A “Worklist” file allows linking an identifier such as a laboratory/sample number to a position of an array within the ArrayStrip. For privacy reasons, arrays should not be identified by patient names. Worklists can be generated using spreadsheet software such as EXCEL (see below) but must be saved in the \*.txt file format that can be imported into the test specific ArrayMate software. Do not use special characters (such as : ; ()[] / \ etc.).

- Create a list with at least three columns that have headers written into the first line. The following headers are obligatory (in this order): position / sample ID / assay ID (Table 1).

- Positions are continuously numbered from 1 to, maximal, 96. Position 1 would correspond to A1, 8 to H1, 9 to A2 and 96 to H12 (Table 2). Do not leave empty lines in the worklist. If you use EXCEL, position numbers should be typed into column A.
- Sample ID are strain/sample/laboratory numbers such as exported from your LIMS (or assigned in any different way). Patient's names should not be used as Sample IDs.
- Assay IDs allow the system to identify the actual test and to correctly use information on layout, spot number and identity etc. *S. aureus* Genotyping Kit has Assay ID: 10248. *This must not be confused with Assay IDs of other tests as this could lead to errors or loss of data.*
- You may add further columns and headers with notes and comments at your convenience. Information from these columns will not appear on the result screens or the Test Report.
- We recommend using a printout of the worklist as template for pipetting.
- Save the worklist as \*.txt file on the memory stick provided together with the ArrayMate.
- To avoid confusion, make sure that worklists are named unambiguously or that worklists from earlier experiments are deleted.



**Table 1:** Example worklist:

Position	sampleID	assayID	comment
1	2013-12345	10248	
2	2013-12346	10248	
3	2013-12347	10248	
4	2013-12348	10248	
5	2013-12349	10248	
6	2013-12350	10248	
7	987654	10248	<i>Isolate referred from Dr. J. Doe.</i>
8	N315	10248	<i>Control strain</i>

**Table 2:** Positions in the 96 well format

	1	2	3	4	5	6	7	8	9	10	11	12
A	1	9	17	25	33	41	49	57	65	73	81	89
B	2	10	18	26	34	42	50	58	66	74	82	90
C	3	11	19	27	35	43	51	59	67	75	83	91
D	4	12	20	28	36	44	52	60	68	76	84	92
E	5	13	21	29	37	45	53	61	69	77	85	93
H	6	14	22	30	38	46	54	62	70	78	86	94
G	7	15	23	31	39	47	55	63	71	79	87	95
H	8	16	24	32	40	48	56	64	72	80	88	96

### Data Acquisition in the ArrayMate Reader

- Insert your memory stick containing the worklist. Use any of the USB ports down to the right side of the ArrayMate.
- Press the button: ; a folder selection dialog will open.
- Select your worklist (path: “My Computer/Removable Disk”).
- Open your selected worklist with “Enter” or the button “Open”.
- Press the button:  (your imported worklist opens in a separate window). Proofread. If the new window is empty or if it was the wrong worklist, repeat the import.
- Press the button “OK”; the worklist window will close.
- Leave the memory stick attached to the ArrayMate if you intend to export *S. aureus* Genotyping Test Reports afterwards.
- Press the button “next” (bottom/right on the screen; reader is opening).
- Carefully insert the appropriate metallic adapter/frame into the ArrayMate. Do not apply any strong force. Assure proper fit, otherwise the images may be out of focus.

- Carefully insert the white frame with the array strips into the metallic adapter. Assure the correct orientation (Position A1 in the frame next to the data matrix barcode on the adapter) and proper fit, otherwise the images may be out of focus.



ArrayStrip frame with inserted strips. Strips are inserted in accordance to the worklist.

*ArrayStrips must be clean. They should not contain any liquids by now. Barcodes must be clean. There must be no lids on the wells that are to be analysed (however, unused wells should remain capped).*

- Press the button “Next” (bottom/right on the screen; reader is closing, analysis program starts, it takes ca. 2-10 min dependent on the number of strips; reader takes images and automatically analyses the data). The progress of the reading is indicated by the following symbols:

photographed:



in analysis:



ready:



- The reader indicates the end of the entire process with an acoustic signal (beep).
- Press the button “Next” (bottom/right on the screen; reader is opening).
- Remove the white frame with the ArrayStrip(s).
- Press the button “Next” (bottom/right on the screen; reader is closing).




## Results

On the left hand site of the screen you will see a list showing all runs stored on the ArrayMate's hard disk. A run contains the results from all arrays analysed together within one frame. If this list is not visible:


- Press the button "Archive" (left hand) and activate the Flag "Browse" (top left).

The runs are organised like folders in "Windows Explorer" and, by default, named according to the date of acquisition.


Example: there is one reading in this archive:

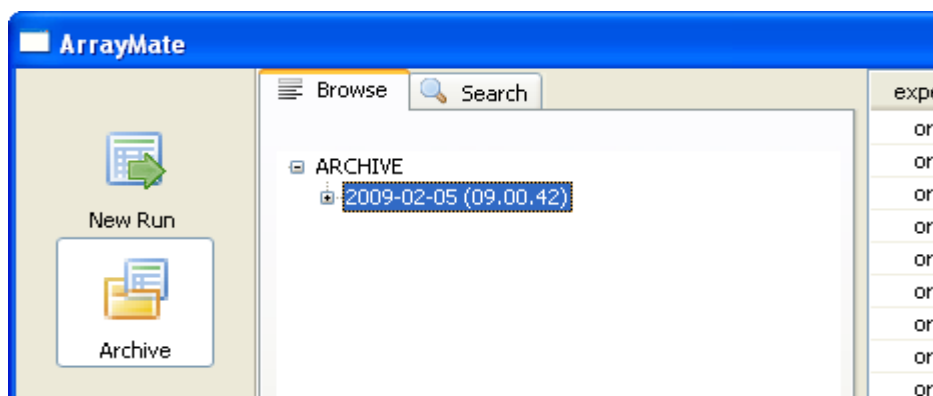

If you click on the plus symbol left on the run name, the folder opens and you will see a list of the individual arrays alphabetically ordered by Sample ID.

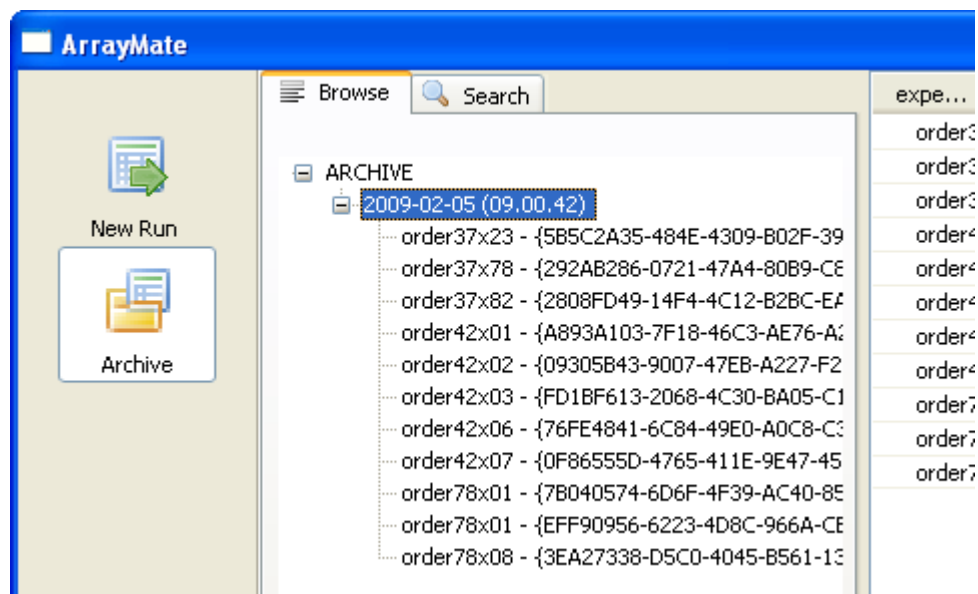

Click on a Sample ID and the *S. aureus* Genotyping Test Report for this array is shown in the window on the right:

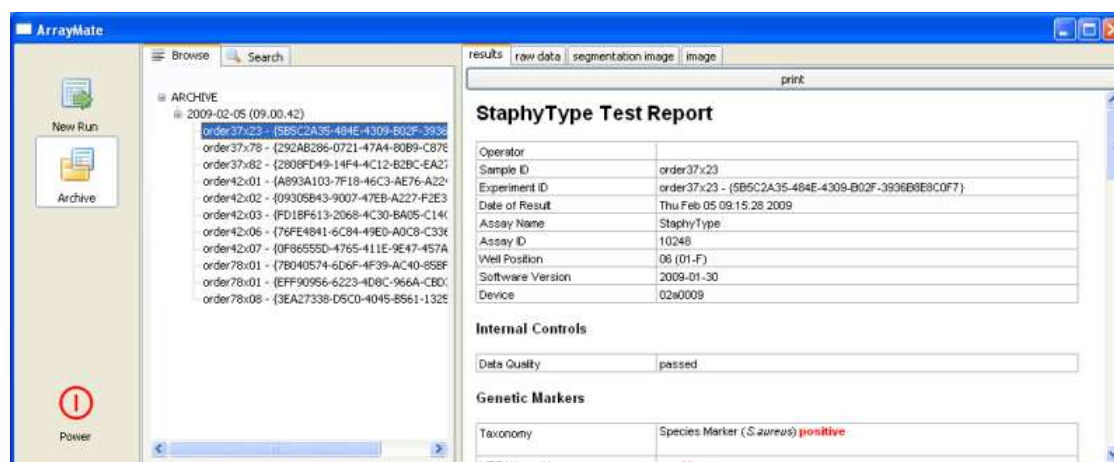

## Export of *S. aureus* Genotyping Test Reports

Two result files in html format will be generated. The shorter one will include a summary on typing information.

This includes the clonal complex affiliation as derived from the overall hybridisation pattern and the strain affiliation as defined by clonal complex affiliation, presence/absence and type of SCCmec elements and presence/absence of PVL or other relevant markers.

MLST sequence types and *spa* types known to be associated with this strain are also displayed. Note that this information is derived from a database search (see also Appendix 3), not from an actual experiment. Furthermore, results for virulence markers and genes associated to antibiotic resistance are listed.

A longer (“result\_B.res.html”) html result sheets will show information on all probes.


Possible error messages that might occur in these reports will be explained below (see Troubleshooting).

Other files that are generated and that can be exported include

- a \*.txt file with the raw measurements,
- an image file (\*.bmp) with the actual photo of the array,
- a second image file (\*.png) in which the coordinate grid is superimposed and the recognised spots are circled and
- \*.xml files that contains the same information as the html result sheets for future export into databases etc.

Note: only complete runs can be exported. The export of individual *S. aureus* Genotyping Test Reports is not possible

- Right-Click on the reading (a menu appears with the option “Export Run Reports”).
- Right-Click on “Export Run Reports” (a file browser opens).

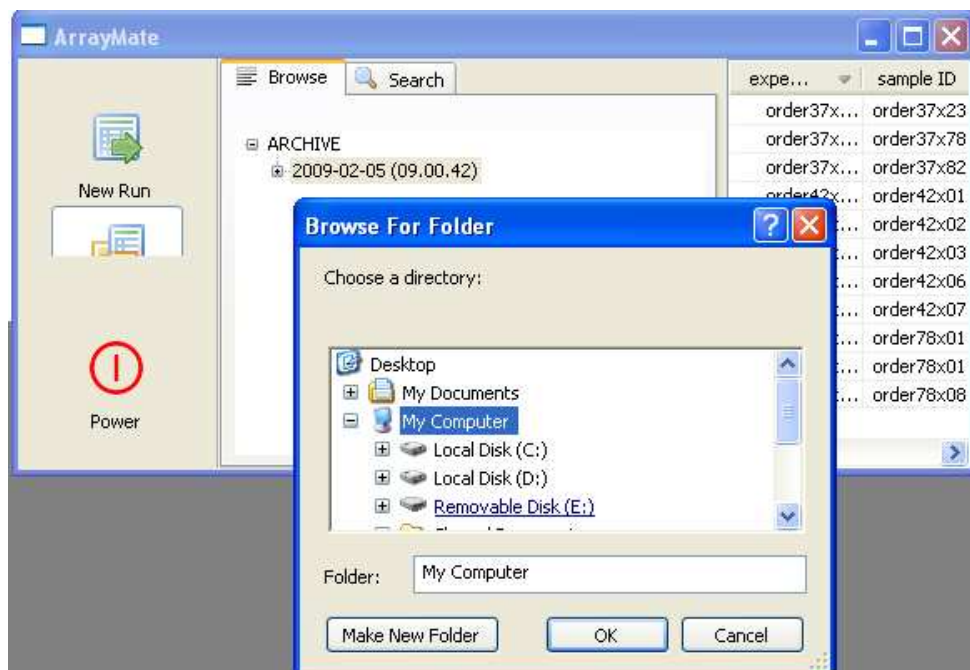

- Click on “My Computer”, then on “Removable Disk” and choose the folder where to store or click on the button “Make New Folder” (on the bottom; a new folder icon appears).
- Rename the new folder (e.g. with the experiment name or date).
- Click on the “Ok” button (data are exported now into the new folder on your memory stick).
- Do NOT remove the memory stick as long as the hourglass symbol is visible.
- Switch off the device by clicking on the “Power”-button (left/down on the screen): ⓘ
- Switch off the Screen. There is no need to physically switch off the ArrayMate.

## **TROUBLESHOOTING**

In case of trouble always make sure that reagents are within the recommended shelf-life and stored in the appropriate way.

In case of trouble we are always happy to support. Please contact [cct.home@clondiag.com](mailto:cct.home@clondiag.com) (or +49 (0) 3641 3111 0), and please include a description of the problem as well as the array images (\*.bmp files) in question.

Please also see Appendix 2 for sample images.

### **Staining Control**

A staining control is included to check whether possible problems originate from the hybridisation or the staining procedure. If the staining control has “Failed” proceed as follows:

Horseradish peroxidase conjugate may have degraded during storage. Add 1 µL buffer C3/C4 to 9 µL D1 (substrate). If the solution turns green within 3-5 seconds, the horseradish peroxidase still has sufficient enzymatic activity.

Enzymatic reaction is inhibited by carryover of buffer C1. Ensure proper washing of the wells to remove all of Buffer C1 prior to adding horseradish peroxidase conjugate.

If the Staining Control has “Passed”, refer to the following hints.

### **Image Quality**

In case of poor image quality we recommend to re-check DNA quantity and quality first by loading leftover DNA on an agarose gel.

In order to determine whether any problems originated from the DNA preparation, perform an experiment with the CM / Control material. This is DNA from the reference strain N315 (GenBank accession number BA000018) and should be identified by the assay as “ST5-MRSA-II [tst1+], New York-Japan Clone”. If the control experiment yields a valid result and a correct identification, there was probably an issue with DNA preparation. If the control experiment also fails, an error affecting later steps or a degradation of reagents from later steps is likely.

### **DNA Quality**

The amount of DNA is crucial because of the linear kinetics of amplification (see Introduction). DNA should be free of RNA, as free RNA reduces the efficiency of amplification and labelling by effectively removing primer from the reaction mix due to competitive hybridisation. A260 readings will cover RNA and other contaminants as well. Therefore pure DNA preparations without RNA contaminations are a prerequisite for proper DNA concentration measurement. RNase treatment prior to A260 reading therefore is necessary (component A2 contains RNase).

DNA must be unfragmented, as fragmentation reduces the efficiency of amplification and labelling due to the distance between primer and probe binding sites. DNA should for this reason not be prepared by disrupting *S. aureus* cells using bead beaters, ultrasonication or aggressive chemicals such as in alkaline lysis protocols. We made good experiences with the manual QIAGEN DNeasy kit and the automated device EZ1.

DNA must be free of any traces of ethanol, as ethanol strongly influences the amplification. It is possible to heat the sample prior to adding it to the labelling mix (5-10 minutes at 70°C). Some problems with samples from the Qiagen EZ1 device for example were resolved after heating the samples (see above).

### **Physical Damage to the Array**

Scratching of the array surface with a pipette tip can lead to the damage of array spots that prohibits the acquisition of a valid signal. In this case the respective marker is not assigned as “negative”, but instead the message “none” appears next to the marker name.

### **Ambiguous Results**

Besides a “positive” or “negative” result for the individual markers on the *S. aureus* Genotyping Test Report, the result can also be “ambiguous”.

In cases affecting resistance genes or virulence factors, no definitive answer with regard to this specific marker can be given. This can be caused by poor sample quality, poor signal quality

and, especially in some resistance-associated genes such *aacA-aphD*, by the presence of plasmids in low copy numbers.

Please note, that for some markers, for which allelic variants were to be discriminated (*bbp*, *clfA*, *clfB* and *fnbB* as well as some *set/ssl* genes, *isaB*, *mprF* and *isdA*), a different approach for analysis was used than for resistance genes or virulence factors. In these genes, alleles that differ only in single nucleotides are recognised. For the sake of the creation of identifiable clonal complex-specific patterns, only the probe with the strongest signal value is regarded as positive, provided that it exceeded breakpoint. All other allele-specific probes are then regarded as ambiguous or, if below the breakpoint, as negative. Therefore it is perfectly normal, if a number of allele-specific probes for these genes yield “ambiguous” signals. The presence or absence of these genes is indicated by fields such as, e.g., “*clfA* (total)” which are summaries for all probes related to the respective gene.

### **Report Unavailable**

If the ArrayMate indicates that no report is available for an array (or multiple arrays on one strip), please check that the strip was positioned properly into the frame. Scratches or drops of condensed water might render the barcode identifier unreadable, please wipe it carefully or try to manually identify the test.

If no obvious reason for the fault can be discovered, please contact the technical service.

### **Error Messages in Result Sheets**

Please compare Appendix 2 for images. If strains cannot be identified, error messages are displayed instead of the short html result sheet. In order to facilitate searching for the cause of the error, the long (“result\_B.res.html”) html result sheet will be generated although it might be faulty. However, it might give a hint what the cause of the problem was.

One possible error message is: *“Identification is not possible. This could be due to technical issues such as poor signal quality, overstaining or to contamination. Please re-clone the culture, confirm its identity as Staphylococcus aureus and its purity, and repeat the experiment.”*

*Identification is also not possible for strains that are not covered by the database. If this is likely (i.e., if your isolate is repeatedly un-identifiable or if you have additional typing data suggesting an unknown strain), please submit the array image and/or the isolate in question to Alere Technologies.”* This will appear for instance when the pattern is entirely irregular or if mutually exclusive alleles are detected simultaneously. The long (“result\_B.res.html”) html result sheets might show in the latter case that several *agr* types or capsule types 5 and 8 were detected in one sample. This can be caused by massive unspecific staining or by contamination / mixed culture. Re-clone and repeat. If this message was prompted by a technically faultless experiment, and if contamination can be ruled out by repeated cloning, please submit the picture and/or the strain for further analysis. It might be an unknown strain that cannot be identified because it was not included into the database. If this is the case we will use multilocus sequence typing (MLST) for further characterisation and might include this strain into future database updates.

This error message in the short result sheet accompanied by positive signals *only* for resistance and SCCmec associated genes indicates the presence of another staphylococcal species (*Staph. epidermidis*, *Staph. haemolyticus* etc.). The long (“result\_B.res.html”) html result sheet should provide this information, occasionally a faulty identification as “*Staph. argenteus*” lineage (CC75), albeit at a low score might, occur.

Another error message “*An assignment to a strain is not possible, although the clonal complex was recognised. This might be caused by technical issues such as poor signal quality, overstaining or contamination. The isolate could also represent a new strain within a known clonal complex, i.e., a strain carrying an unusual SCCmec element or an unusual set of virulence genes. If this appears to be the case, please submit the array image and/or the isolate in question to Alere Technologies*” might appear instead of the typing information in an otherwise normal result file. This could indicate an unusual SCCmec element or an unusual presence of virulence genes, such as of PVL in a lineage where it has not been observed before. A contamination, e.g., by SCC-bearing coagulase-negatives needs to be ruled out. Re-clone and repeat. If this message was prompted by a technically faultless experiment, and if contamination can be ruled out by repeated cloning, please submit the picture and/or the strain for further analysis.



## **ADDITIONAL INFORMATION**

### **Warranty**

Alere guarantees the performance as described in this manual. Usage of the Kit was successfully tested at ambient temperatures up to 37°C, a guarantee is limited to ambient temperatures in the laboratory between 18°C and 28°C. Kit components comprise the arrays and their caps, the Lysis Enhancer, the reagents for DNA labelling and for detection of labelled DNA products on the array, the ArrayMate reader and its software. In case one of these components fails within the expiry date due to other reason than misuse, contact Alere for replacement or refund. Terms and conditions apply.

If you have any problem or question, please contact the technical service.

### **Disclaimer**

***This system is for research use only.***

We do not accept any liability for damages caused by misuse. Misuse comprises, especially but not exclusively, of a use of the system for the detection of resistance genes in order to predict phenotypic antibiotic resistances or susceptibilities for the guidance of an antibiotic chemotherapy.

*Since resistances might be caused by genes or mutations not covered by this array or by hitherto unknown genes or mutations, any antibiotic chemotherapy MUST be guided by phenotypic susceptibility tests.*

Furthermore, we do not accept any liability for damages caused by inappropriate use of the device as a personal computer, for instance related to the use of additional software, to network connections, or to a breach of privacy related to the storage of confidential information (such as names of patients from whom *S. aureus* was isolated) on its hard disk and/or to the use of external storage devices that might be contaminated with spyware.

### **Quality Control**

Each batch is stringently tested with the use of standard *S. aureus* DNA preparations for good performance and correctness of results.

### **List of Components for Separate Order**

If required, these reagents for the *S. aureus* Genotyping Kit may be ordered separately:

Component	Name	Amount	Cat#	Storage
A1	Lysis Buffer	30 ml	245101000	18-28 °C
A2	Lysis Enhancer	96 units	245102000	18-28 °C
B1 <sup>ST</sup>	Labelling Buffer/Master Mix	700 µl	245103000	2-8 °C
B2	Labelling Enzyme	20 µl	245104000	2-8 °C
C1	Hybridisation Buffer	30 ml	245105000	18-28 °C
C2	Washing Buffer 1	120 ml	245106000	18-28 °C
C3	HRP Conjugate 100x	200 µl	245107000	2-8 °C
C4	Conjugate Buffer	30 ml	245108000	18-28 °C
C5	Washing Buffer 2	120 ml	245109000	18-28 °C
D1	HRP Substrate	15 ml	245110000	2-8 °C
ArrayStrips	ArrayStrips for <i>S. aureus</i> 2.0	12 St	240008560	15-28 °C
StripCaps	Covers for unused arrays	12 St	245112000	18-28 °C
CM	Control Material (N315 DNA)	30 µl	245111000	2-8 °C

For pricing please contact your local representative or our customer service, respectively.

### **Legal Manufacturer**

Alere Technologies GmbH

Loebstedter Str. 103-105

07749 Jena, Germany

### **Contact**

If you require any further information on this product please contact [cct.home@clondiag.com](mailto:cct.home@clondiag.com)

## **LITERATURE**

Literature quoted in this manual:

- [1] Monecke S, Coombs G, Shore AC, Coleman DC, Akpaka P, Borg M, Chow H, Ip M, Jatzwauk L, Jonas D, Kadlec K, Kearns A, Laurent F, O'Brien FG, Pearson J, Ruppelt A, Schwarz S, Scicluna E, Slickers P, Tan H-L, Weber S, Ehricht R (2011) A Field Guide to Pandemic, Epidemic and Sporadic Clones of Methicillin-Resistant *Staphylococcus aureus*. PLoS One 6 (4):e17936
- [2] Monecke S, Slickers P, Ehricht R (2008) Assignment of *Staphylococcus aureus* isolates to clonal complexes based on microarray analysis and pattern recognition. FEMS Immunol Med Microbiol 53:237–251

For further literature please refer to:

<http://alere-technologies.com/en/science-technologies/publications/saureus.html>

## **UPDATES & SOFTWARE**

Notifications on database/software updates and freeware tools can be found at:

<http://alere-technologies.com/en/science-technologies/publications/downloads.html>.


and/or <http://alere-technologies.com/en/news.html>.

Currently available freeware programs are:

- “spa type mapper” for the analysis of spa sequences
- “Alere S aureus Results Collector” for the conversion of multiple \*result.xml files from the ArrayMate into spreadsheet tables. This should make it easier to compare isolates or to determine relative abundances of genes or strains etc.


## APPENDIX 1 - FLOW CHART

**Quantifoil protocol.** The figure on this page summarises the test procedure for the thermoshaker BioShake iQ by Quantifoil. Please always refer to the text section of this manual for further important details.




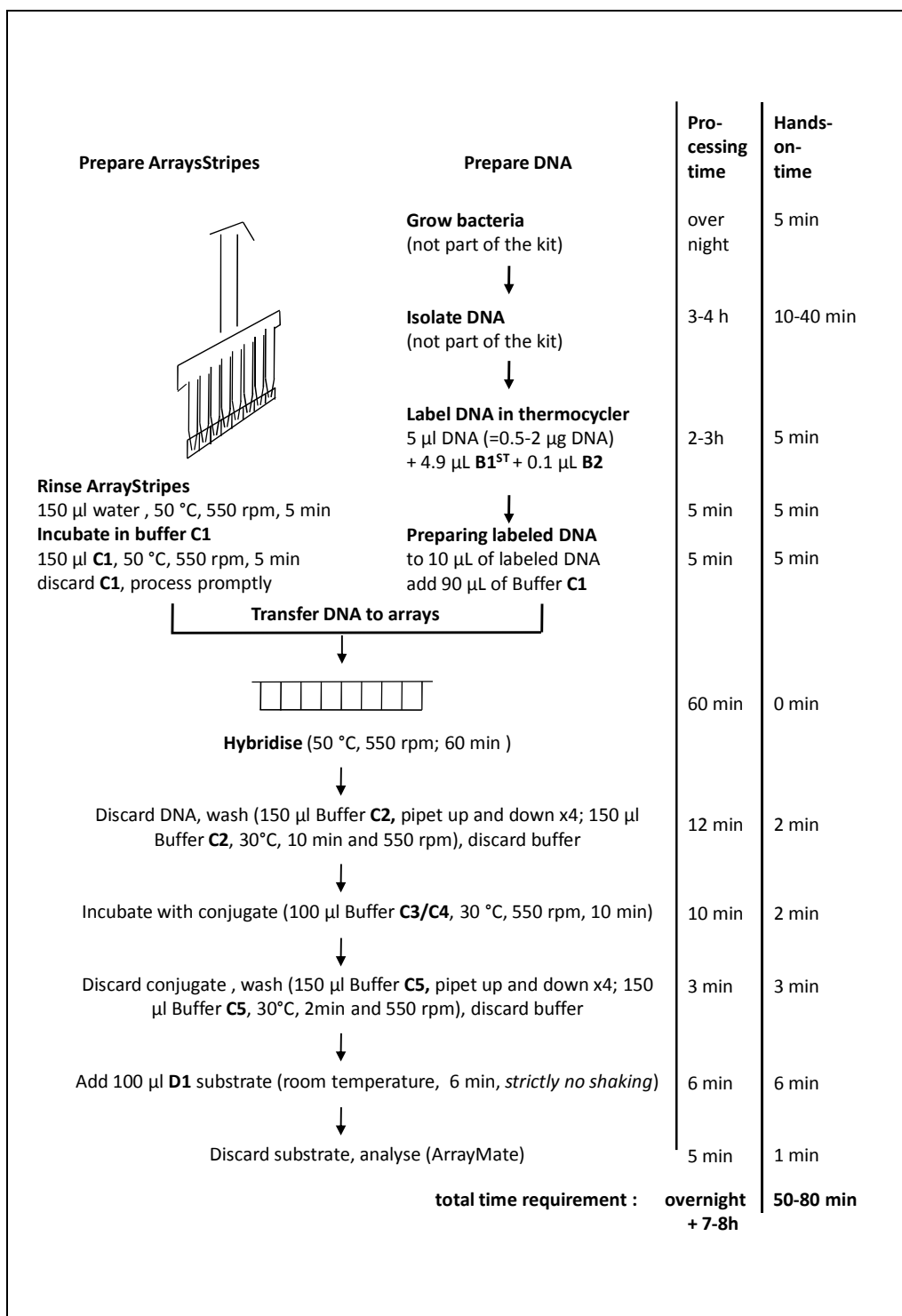


## APPENDIX 1 - FLOW CHART

**Eppendorf protocol.** The figure on this page summarises an adapted test procedure for the thermoshaker Thermomixer Comfort by Eppendorf. Please always refer to the text section of this manual for further important details.

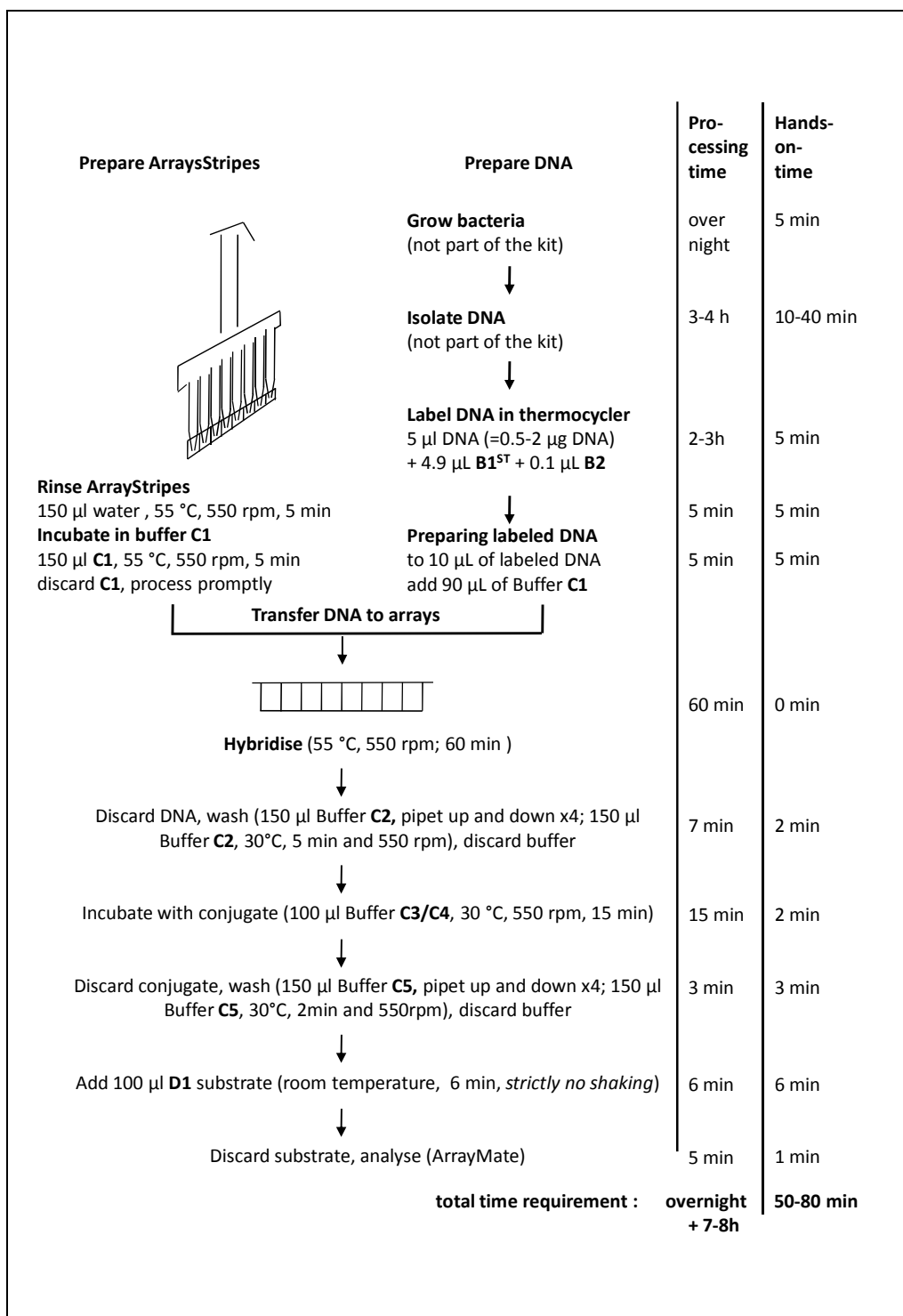

## APPENDIX 2 – IMAGES FOR TROUBLESHOOTING

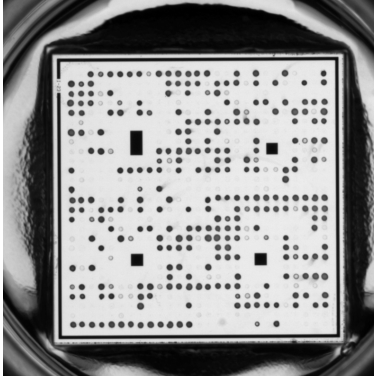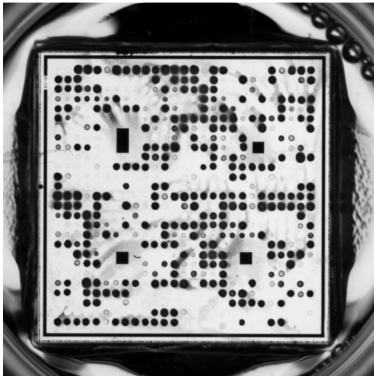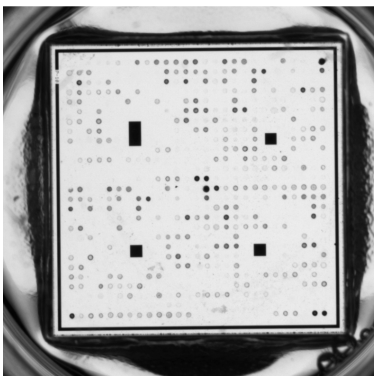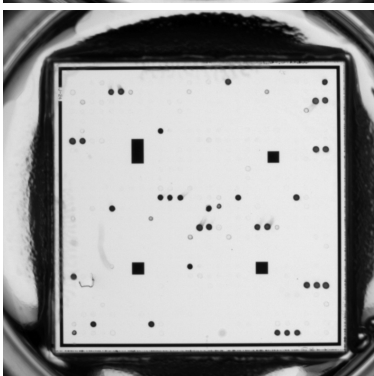Image	Comment	Result sheets:
	A technically faultless, valid experiment.	Valid results, no error messages.
	This image is overstained. The experiment should be repeated.	There might be no error messages although individual probes might yield false-positives. The error message <i>"An assignment to a strain is not possible, although the clonal complex was recognised. This might be caused by technical issues such as poor signal quality, overstaining or contamination. ..."</i> might appear if false-positives hinder strain identification.
	This image is poor. This could be due to low DNA concentration, fragmented DNA, ethanol trace contaminations in DNA sample or expired reagents. The experiment should be repeated with a new DNA preparation. If this also fails, try an experiment with N315 control DNA (CM).	The error message <i>"An assignment to a strain is not possible, although the clonal complex was recognised. This might be caused by technical issues such as poor signal quality, overstaining or contamination. The isolate could also represent a new strain within a known clonal complex ..."</i> might appear in the short file. The long file will yield an approximate identification, but individual probes might yield false-negative results.
	Species other than <i>S. aureus</i> tested. Check identification by other means.	Error message in the short html file. The long file yields <i>"Coagulase-negative Staphylococci, other bacteria, or very poor signal quality. Check identification by biochemical means or MALDI-TOF or repeat experiment"</i> and shows <u>positive results only for resistance genes and/or genes associated with SCCmec</u> .

## APPENDIX 3 – TARGET GENES

SPECIES MARKER	domain 1 of 23S-rRNA	<i>rrnD1 (S. aureus)</i>
	glyceraldehyde 3-phosphate dehydrogenase, locus 1	<i>gapA</i>
	katalase A	<i>katA</i>
	coagulase	<i>coA</i>
	thermostable extracellular nuclease	<i>nuc1</i>
	staphylococcal protein A	<i>spa</i>
	IgG-binding protein	<i>sbi</i>
REGULATORY GENES	staphylococcal accessory regulator A	<i>sarA</i>
	histidine protein kinase, sae locus	<i>saeS</i>
	sensor protein	<i>vraS</i>
	accessory gene regulator allele I	<i>agrI (total)</i>
		<i>agrB-I</i>
		<i>agrC-I</i>
		<i>agrD-I</i>
	accessory gene regulator allele II	<i>agrII (total)</i>
		<i>agrB-II</i>
		<i>agrC-II</i>
		<i>agrD-II</i>
	accessory gene regulator allele III	<i>agrIII (total)</i>
		<i>agrB-III</i>
		<i>agrC-III</i>
		<i>agrD-III</i>
	accessory gene regulator allele IV	<i>agrIV (total)</i>
		<i>agrB-IV</i>
		<i>agrC-IV</i>
	haemolysin delta	<i>hld</i>
METHICILLIN RESISTANCE AND SCCmec TYPING	alternate penicillin binding protein 2, defining MRSA	<i>mecA</i>
	truncated signal transducer protein MecR1	<i>delta_mecR</i>
	glycerophosphoryl diester phosphodiesterase, associated with mecA	<i>ugpQ</i>
	cassette chromosome recombinase genes A-1	<i>ccrA-1</i>
	cassette chromosome recombinase genes B-1	<i>ccrB-1</i>
	plasmin-sensitive surface protein	<i>plsSCC (COL)</i>
	hypothetical protein from SCCmec elements	<i>Q9XB68-dcs</i>
	cassette chromosome recombinase gene A-2	<i>ccrA-2</i>
	cassette chromosome recombinase gene B-2	<i>ccrB-2</i>
	potassium-translocating ATPase A, chain 2	<i>kdpA-SCC</i>
	potassium-transporting ATPase B, chain 1	<i>kdpB-SCC</i>
	potassium-translocating ATPase C, chain 2	<i>kdpC-SCC</i>
	sensor kinase protein	<i>kdpD-SCC</i>
	KDP operon transcriptional regulatory protein	<i>kdpE-SCC</i>
	methicillin-resistance gene regulatory protein	<i>mecl</i>
	signal transducer protein MecR1	<i>mecR</i>
	homolog of xylose repressor, associated with SCCmec-elements	<i>xylR</i>

	cassette chromosome recombinase gene A-3	<i>ccrA-3</i>
	cassette chromosome recombinase gene B-3	<i>ccrB-3</i>
	mercury resistance gene operon, Hg(II) reductase	<i>merA</i>
	mercury resistance gene operon, alkylmercury lyase	<i>merB</i>
	Putative protein, homologue to cassette chromosome recombinase A genes	<i>ccrAA (MRSZ47)_probe 1</i> <i>ccrAA (MRSZ47)_probe 2</i>
	cassette chromosome recombinase gene C	<i>ccrC (85-2082)</i>
	cassette chromosome recombinase gene A-4	<i>ccrA-4</i>
	cassette chromosome recombinase gene B-4	<i>ccrB-4</i>
RESISTANCE : PENICILLINASE	beta-lactamase gene	<i>blaZ</i>
	beta lactamase repressor (inhibitor)	<i>blal</i>
	beta-lactamase regulatory protein	<i>blaR</i>
RESISTANCE : MLS- ANTIBIOTICS	rRNA methyltransferase associated with macrolide/lincosamide resistance	<i>erm(A)</i>
	rRNA methyltransferase associated with macrolide/lincosamide resistance	<i>erm(B)</i>
	rRNA methyltransferase associated with macrolide/lincosamide resistance	<i>erm(C)</i>
	lincosaminide nucleotidyltransferase (=linA)	<i>lnu(A)</i>
	macrolide efflux pump	<i>msr(A)</i>
	macrolide efflux protein A	<i>mef(A)</i>
	macrolide phosphotransferase II (=mpbBM)	<i>mph(C)</i>
	virginiamycin A acetyltransferase	<i>vat(A)</i>
	acetyltransferase inactivating streptogramin A	<i>vat(B)</i>
	ABC transporter conferring resistance to streptogramin A and related compounds	<i>vga(A)</i>
	vga(A) allele from strain BM 3327	<i>vga(A) (BM 3327)</i>
	virginiamycin B hydrolase (=vgb)	<i>vgb(A)</i>
RESISTANCE : AMINOGLYCOSIDES	aminoglycoside adenyl-/phosphotransferase (gentamicin, tobramycin)	<i>aacA-aphD</i>
	aminoglycoside adenyltransferase (neo-/ kanamycin, tobramycin)	<i>aadD</i>
	aminoglycoside phosphotransferase (neo-/ kanamycin)	<i>aphA3</i>
RESISTANCE : MISCELLANEOUS GENES	streptothricin acetyltransferase	<i>sat</i>
	dihydrofolate reductase mediating trimethoprim resistance (=dfrA)	<i>dfrS1</i>
	fusidic acid resistance gene (= far1)	<i>fusB</i>
	fusidic acid resistance gene (= Q6GD50)	<i>fusC</i>
	isoleucyl-tRNA synthetase associated with mupirocin resistance (=mupR)	<i>mupA</i>
	tetracycline efflux protein	<i>tet(K)</i>
	ribosomal protection protein associated with tetracycline resistance	<i>tet(M)</i>
	chloramphenicol acetyltransferase	<i>cat (total)</i> <i>cat (pC221)</i> <i>cat (pC223)</i> <i>cat (pMC524)</i> <i>cat (pSBK203R)</i>
	23S rRNA methyltransferase (phenicols, lincosamides, oxazolidinones, pleuromutilins, streptogramin A )	<i>cfr</i>



	chloramphenicol/florfenicol exporter	<i>fexA</i>
	metallothiol transferase	<i>fosB</i>
		<i>fosB (plasmid)</i>
RESISTANCE : EFFLUX SYSTEMS	quaternary ammonium compound / multidrug efflux protein C	<i>qacA</i>
	quaternary ammonium compound / multidrug efflux protein A	<i>qacC (total)</i>
		<i>qacC (consensus)</i>
		<i>qacC (equine)</i>
		<i>qacC (SA5)</i>
		<i>qacC (Ssap)</i>
		<i>qacC (ST94)</i>
	putative transport protein (=tetEfflux)	<i>sdrM</i>
RESISTANCE : GLYCOPEPTIDES	vancomycin resistance gene	<i>vanA</i>
	vancomycin resistance gene from enterococci and Clostridium	<i>vanB</i>
	teicoplanin resistance gene from enterococci	<i>vanZ</i>
VIRULENCE : TOXIC SCHOCK TOXIN	toxic shock syndrome toxin 1	<i>tst1 (consensus)</i>
		<i>tst1 ("human" allele)</i>
		<i>tst1 ("bovine" allele, from RF122)</i>
VIRULENCE : ENTEROTOXINS	enterotoxin A (=entA)	<i>sea</i>
	enterotoxin A, allele from strain 320E	<i>sea (320E)</i>
	enterotoxin A, allele from strain N315 = enterotoxin P	<i>sea (N315)</i>
	enterotoxin B (=entB)	<i>seb</i>
	enterotoxin C (=entC)	<i>sec</i>
	enterotoxin D (=entD)	<i>sed</i>
	enterotoxin E (=entE)	<i>see</i>
	enterotoxin G (=entG)	<i>seg</i>
	enterotoxin H (=entH)	<i>seh</i>
	enterotoxin I (=entI)	<i>sei</i>
	enterotoxin J (=entJ)	<i>sej</i>
	enterotoxin K (=entK)	<i>sek</i>
	enterotoxin L (=entL)	<i>sel</i>
	enterotoxin-like gene/protein M (=sem, entM)	<i>selm</i>
	enterotoxin-like gene/protein N (=sen, entN)	<i>seln (consensus)</i>
		<i>seln (other than RF122)</i>
	enterotoxin-like gene/protein O (=seo, entO)	<i>selo</i>
	enterotoxin gene cluster (seg/i/selm/n/o/u)	<i>egc</i>
	enterotoxin Q (=entQ)	<i>seq</i>
	enterotoxin R (=entR)	<i>ser</i>
	enterotoxin-like gene/protein U (=seu, entU)	<i>selu</i>
	enterotoxin-like protein ORF CM14	<i>ORF CM14_ probe1</i>
	enterotoxin-like protein ORF CM14	<i>ORF CM14_ probe2</i>
VIRULENCE : HLG AND LEUKOCIDINS	haemolysin gamma / leukocidin, component B (F)	<i>lukF</i>
	haemolysin gamma / leukocidin, component C (S)	<i>lukS</i>
	haemolysin gamma / leukocidin, component C (S), allele from ST22 and ST45	<i>lukS (ST22+ST45)</i>
	haemolysin gamma, component A	<i>hlgA</i>

	Panton Valentine leukocidin F component	<i>lukF-PV</i>
	Panton Valentine leukocidin S component	<i>lukS-PV</i>
	F component of leukocidin from ruminants	<i>lukF-PV (P83)</i>
	S component of leukocidin from ruminants	<i>lukM</i>
	leukocidin D component	<i>lukD</i>
	leukocidin E component	<i>lukE</i>
	leukocidin/ haemolysin toxin family protein	<i>lukX</i>
	leukocidin/haemolysin toxin family protein, allele from ST30 and ST45	<i>lukY</i>
	leukocidin/haemolysin toxin family protein	<i>lukY (ST30+ST45)</i>
VIRULENCE : HAEMOLYSINS	putative membrane protein	<i>hl</i>
	haemolysin alpha	<i>hla</i>
	putative membrane protein	<i>hlIII (consensus)</i>
		<i>hlIII (other than RF122)</i>
	haemolysin beta	<i>hly_probe 1</i>
	haemolysin beta	<i>hly_probe 2</i>
	haemolysin beta	<i>hly_probe 3</i>
	haemolysin beta without phage insertion	<i>un-disrupted hly</i>
VIRULENCE : HLB-CONV PHAGES	staphylokinase	<i>sak</i>
	chemotaxis-inhibiting protein (CHIPS)	<i>chp</i>
	staphylococcal complement inhibitor	<i>scn</i>
VIRULENCE : EXFOLIATIVE TOXINS	exfoliative toxin serotype A	<i>etA</i>
	exfoliative toxin serotype B	<i>etB</i>
	exfoliative toxin D	<i>etD</i>
VIRULENCE : EPITHEL. DIFF. INHIB	epidermal cell differentiation inhibitor	<i>edinA</i>
	epidermal cell differentiation inhibitor B	<i>edinB</i>
	epidermal cell differentiation inhibitor C	<i>edinC</i>
VIRULENCE : ACME LOCUS	Arginine Catabolic Mobile Element	<i>ACME cluster</i>
	ACME-locus: arginine deiminase	<i>arcA-SCC</i>
	ACME-locus: ornithincarbamoyltransferase	<i>arcB-SCC</i>
	ACME-locus: carbamyltransferase	<i>arcC-SCC</i>
	ACME-locus: arginine/ornithine-antiporter	<i>arcD-SCC</i>
VIRULENCE : PROTEASES	aureolysin	<i>aur (consensus)</i>
		<i>aur (other than MRSA252)</i>
		<i>aur (MRSA252)</i>
	serinprotease A	<i>splA</i>
	serinprotease B	<i>splB</i>
	serinprotease E	<i>splE</i>
	glutamylendopeptidase	<i>sspA</i>
	staphopain B, protease	<i>sspB</i>
	staphopain A (staphylopain A), protease	<i>sspP (consensus)</i>
		<i>sspP (other than ST93)</i>
VIRULENCE : STAPHYLOCOCCAL SUPERANTIGEN/ ENTEROTOXIN-LIKE GENES (SET/SSL)	staphylococcal exotoxin-like protein/SAg gene homolog, SAUSA300_0370	<i>setC/setx</i>
	staphylococcal superantigen-like protein 1 (probes)	<i>ssl01/set6_probe1_11</i>
		<i>ssl01/set6_probe2_11</i>
		<i>ssl01/set6_probe1_12</i>
		<i>ssl01/set6_probe2_12</i>
		<i>ssl01/set6_probe4_11</i>

	<i>ssl01/set6_probeRF122</i>
staphylococcal superantigen-like protein 1 (interpretation/alleles)	<i>ssl01/set6 (COL)</i>
	<i>ssl01/set6 (Mu50+N315)</i>
	<i>ssl01/set6 (MW2+MSSA476)</i>
	<i>ssl01/set6 (MRSA252)</i>
	<i>ssl01/set6 (RF122)</i>
	<i>ssl01/set6 (other alleles)</i>
staphylococcal superantigen-like protein 2	<i>ssl02/set7</i>
	<i>ssl02/set7 (MRSA252)</i>
staphylococcal superantigen-like protein 3	<i>ssl03/set8_probe 1</i>
	<i>ssl03/set8_probe 2</i>
	<i>ssl03/set8 (MRSA252, SAR0424)</i>
staphylococcal superantigen-like protein 4	<i>ssl04/set9</i>
	<i>ssl04/set9 (MRSA252, SAR0425)</i>
staphylococcal superantigen-like protein 5	<i>ssl05/set3_probe 1</i>
	<i>ssl05/set3 (RF122, probe-611)</i>
	<i>ssl05/set3_probe 2 (612)</i>
	<i>ssl05/set3 (MRSA252)</i>
staphylococcal superantigen-like protein 6	<i>ssl06/set21</i>
	<i>ssl06 (NCTC8325+MW2)</i>
staphylococcal superantigen-like protein 7	<i>ssl07/set1</i>
	<i>ssl07/set1 (MRSA252)</i>
	<i>ssl07/set1 (AF188836)</i>
staphylococcal superantigen-like protein 8	<i>ssl08/set12_probe 1</i>
	<i>ssl08/set12_probe 2</i>
staphylococcal superantigen-like protein 9	<i>ssl09/set5_probe 1</i>
	<i>ssl09/set5_probe 2</i>
	<i>ssl09/set5 (MRSA252)</i>
staphylococcal superantigen-like protein 10	<i>ssl10/set4</i>
	<i>ssl10 (RF122)</i>
	<i>ssl10/set4 (MRSA252)</i>
staphylococcal superantigen-like protein 11	<i>ssl11/set2 (COL)</i>
	<i>ssl11+set2(Mu50+N315)</i>
	<i>ssl11+set2(MW2+MSSA476)</i>
	<i>ssl11/set2 (MRSA252)</i>
staphylococcal exotoxin-like protein, second locus	<i>setB3</i>
	<i>setB3 (MRSA252)</i>
	<i>setB2</i>
	<i>setB2 (MRSA252)</i>
	<i>setB1</i>
CAPSULE- AND BIOFILM- ASSOCIATED GENES	capsule type 1
	<i>cap 1 (total)</i>
	capsular polysaccharide synthesis enzyme
	<i>capH1</i>
	O-antigen polymerase
	<i>capJ1</i>
	capsular polysaccharide biosynthesis protein
	<i>capK1</i>
	capsule type 5
	<i>cap 5 (total)</i>

	capsular polysaccharide synthesis enzyme	<i>capH5</i>
	O-antigen polymerase	<i>capJ5</i>
	capsular polysaccharide biosynthesis protein	<i>capK5</i>
	capsule type 8	<i>cap 8 (total)</i>
	capsular polysaccharide synthesis enzyme	<i>capH8</i>
	capsular polysaccharide biosynthesis protein	<i>capI8</i>
	O-antigen polymerase	<i>capJ8</i>
	capsular polysaccharide biosynthesis protein	<i>capK8</i>
	intercellular adhesion protein A	<i>icaA</i>
	intercellular adhesion protein C	<i>icaC</i>
	biofilm PIA synthesis protein D	<i>icaD</i>
	surface protein involved in biofilm formation	<i>bap</i>
ADHAESION FACTORS / GENES ENCODING MICROBIAL SURFACE COMPONENTS RECOGNIZING ADHESIVE MATRIX MOLECULES (MSCRAMM GENES)	bone sialoprotein-binding protein	<i>bbp (total)</i>
		<i>bbp (consensus)</i>
		<i>bbp (COL+MW2)</i>
		<i>bbp (MRSA252)</i>
		<i>bbp (Mu50)</i>
		<i>bbp (RF122)</i>
	clumping factor A	<i>bbp (ST45)</i>
		<i>clfA (total)</i>
		<i>clfA (consensus)</i>
		<i>clfA (COL+RF122)</i>
		<i>clfA (MRSA252)</i>
		<i>clfA (Mu50+MW2)</i>
	clumping factor B	<i>clfB (total)</i>
		<i>clfB (consensus)</i>
		<i>clfB (COL+Mu50)</i>
		<i>clfB (MW2)</i>
		<i>clfB (RF122)</i>
	collagen-binding adhesin	<i>cna</i>
	cell wall associated fibronectin-binding protein	<i>ebh (consensus)</i>
	cell surface elastin binding protein	<i>ebpS (total)</i>
		<i>ebpS_probe 612</i>
		<i>ebpS_probe 614</i>
		<i>ebpS (01-1111)</i>
		<i>ebpS (COL)</i>
	enolase	<i>eno</i>
	fibrinogen binding protein (19 kDa)	<i>fib</i>
		<i>fib (MRSA252)</i>
	fibronectin-binding protein A	<i>fnbA (total)</i>
		<i>fnbA (consensus)</i>
		<i>fnbA (COL)</i>
		<i>fnbA (MRSA252)</i>
		<i>fnbA (Mu50+MW2)</i>
	fibronectin-binding protein B	<i>fnbA (RF122)</i>
		<i>fnbB (total)</i>
		<i>fnbB (COL)</i>
		<i>fnbB (COL+Mu50+MW2)</i>

		<i>fnbB</i> (Mu50)
		<i>fnbB</i> (MW2)
		<i>fnbB</i> (ST15)
		<i>fnbB</i> (ST45-2)
	major histocompatibility complex class II analog protein (=Extracellular adherence protein, eap)	<i>map</i> (total)
		<i>map</i> (COL)
		<i>map</i> (MRSA252)
		<i>map</i> (Mu50+MW2)
	Staphylococcus aureus surface protein G	<i>sasG</i> (total)
		<i>sasG</i> (COL+Mu50)
		<i>sasG</i> (MW2)
		<i>sasG</i> (other than MRSA252+RF122)
	Ser-Asp rich fibrinogen-/ bone sialoprotein-binding protein C	<i>sdrC</i> (total)
		<i>sdrC</i> (consensus)
		<i>sdrC</i> (B1)
		<i>sdrC</i> (COL)
		<i>sdrC</i> (Mu50)
		<i>sdrC</i> (MW2+MRSA252+RF122)
		<i>sdrC</i> (other than MRSA252+RF122)
	Ser-Asp rich fibrinogen-/ bone sialoprotein-binding protein D	<i>sdrD</i> (total)
		<i>sdrD</i> (consensus)
		<i>sdrD</i> (COL+MW2)
		<i>sdrD</i> (Mu50)
		<i>sdrD</i> (other)
	van Willebrand factor binding protein	<i>vwb</i> (total)
		<i>vwb</i> (consensus)
		<i>vwb</i> (COL+MW2)
		<i>vwb</i> (MRSA252)
		<i>vwb</i> (Mu50)
		<i>vwb</i> (RF122)
IMMUNODOMINANT ANTIGEN B	immunodominant antigen B	<i>isaB</i>
		<i>isaB</i> (MRSA252)
DEFENSIN RESISTANCE	defensin resistance gene protein	<i>mprF</i> (COL+MW2)
		<i>mprF</i> (Mu50+MRSA252)
TRANSFERRIN BINDING PROT	transferrin-binding protein	<i>isdA</i> (consensus)
		<i>isdA</i> (MRSA252)
		<i>isdA</i> (other than MRSA252 )
PUTATIVE TRANSPORTER	hypothetical protein, similar to integral membrane protein LmrP	<i>lmrP</i> (other than RF122)_probe1
		<i>lmrP</i> (other than RF122)_probe2
		<i>lmrP</i> (RF122)_probe1
		<i>lmrP</i> (RF122)_probe2
TYPE I RESTRICTION- MODIFICATION SYSTEM, SINGLE SEQUENCE SPECIFICITY PROTEIN	type I site-specific deoxyribonuclease subunit, 1st locus	<i>hdsS1</i> (RF122)
	type I site-specific deoxyribonuclease subunit, 2nd locus	<i>hdsS2</i> (Mu50+N315+COL+USA300+ NCTC8325)

		<i>hsdS2</i> (MW2+MSSA476)
		<i>hsdS2</i> (RF122)
		<i>hsdS2</i> (MRSA252)
	type I site-specific deoxyribonuclease subunit, 3rd locus	<i>hsdS3</i> (all other than RF122+MRSA252)
		<i>hsdS3</i> (COL+USA300+NCTC8325+MW2+MSSA476+RF122)
		<i>hsdS3</i> (Mu50+N315)
		<i>hsdS3</i> (CC51+ MRSA252)
		<i>hsdS3</i> (MRSA252)
	type I site-specific deoxyribonuclease subunit, unknown locus	<i>hsdSx</i> (CC25)
		<i>hsdSx</i> (CC15)
		<i>hsdSx</i> (etd)
MISCELLANEOUS GENES	hypothetical protein, located next to serine protease operon	<i>Q2FXC0</i>
	unspecific efflux/transporter	<i>Q2YUB3</i>
	hypothetical protein	<i>Q7A4X2</i>
HYALURONATE LYASE	hyaluronate lyase, first / second locus	<i>hysA1</i> (MRSA252)
		<i>hysA1</i> (MRSA252+RF122) and/or <i>hysA2</i> (consensus)
		<i>hysA1</i> (MRSA252+RF122) and/or <i>hysA2</i> (COL+USA300)
	hyaluronate lyase, second locus	<i>hysA2</i> (all other than MRSA252)
		<i>hysA2</i> (COL+USA300+NCTC8325)
		<i>hysA2</i> (all other than COL+USA300+NCTC8325)_ <i>probe1</i>
		<i>hysA2</i> (all other than COL+USA300+NCTC8325)_ <i>probe1</i>
		<i>hysA2</i> (MRSA252)

NOTE: “(consensus)” indicates true consensus probes while “(total)” indicates a summary for all probes for a given gene to show on one glance whether this gene is present, in any known allele, or not.

## **APPENDIX 4 – TYPING INFORMATION**

### **Definitions & Explanations**

The displayed result will yield following typing information:

- Strain assignment, as determined by overall profile and by preset definitions for strains. Strains are always defined by clonal complex affiliation (see below), absence or presence of *mecA* and by *SCCmec* type as well as absence or presence of PVL. Widely or historically recognised strains might also be defined based on the absence or presence of additional characteristic genes.
- MLST clonal complex affiliation, as determined by overall profile and by preset definitions including capsule types and *agr* groups. Isolates can be assigned to clonal complexes as defined by multilocus sequence typing. Analysis of hybridisation patterns cannot discriminate sequence types which differ only in single point mutations affecting MLST genes (e.g., ST5 and ST225, or ST59 and ST952). However, there are also sequence types which originate from chromosomal replacements. Examples are CC8/ST239 or CC30/ST34. As these events result in different hybridisation patterns, such STs can be easily identified. Some other STs are also clearly different from parental CCs although recombination is not (yet?) proven. In such cases, ST affiliation might also be displayed (for instance for CC9/ST834)
- Strain synonyms. These are listed if they can be unambiguously attached to strains as defined above. If you use local designations for strains that you want to be included, please contact [stefan.monecke@clondiag.com](mailto:stefan.monecke@clondiag.com)
- Assignment score. This is a score for the similarity to the average hybridisation result for a given strain/CC. Scores below 88% exclude reliable strain identification, and could be attributed either to technical reasons or to the presence of a yet unknown strain. The “short” html file will display an error message; the “long” html file will display the most similar strain although that identification might be faulty. A value of 100% is unlikely because of the mobility of many genes in *S. aureus*.

## List of Currently Recognised Strains

This list corresponds to the software version as of January 2012. If you have array images of a strain not yet covered, please contact [stefan.monecke@clondiag.com](mailto:stefan.monecke@clondiag.com).

Clonal complex	Strain
CC1	CC1-MSSA
	CC1-MSSA [PVL+]
	CC1-MSSA-SCCfus
	CC1-MSSA-SCCfus [PVL+]
	CC1-MRSA-IV [PVL+], USA400
	CC1-MRSA-IV&SCCfus [PVL+]
	CC1-MRSA-IV&SCCfus, WA MRSA-1/45
	CC1-MRSA-IV, WA MRSA-1/57
	CC1-MRSA-V
	CC1-MRSA-V [PVL+]
	CC1-MRSA-V&SCCfus
	CC1-MRSA-V&SCCfus [PVL+]
CC1 (ST567)	ST567-MSSA [PVL+]
CC1 (ST573/772)	ST573/772-MSSA
	ST573/772-MSSA [PVL+]
	ST573-MRSA-V, WA MRSA-10
	ST772-MRSA-V
	ST772-MRSA-V [PVL+], Bengal Bay Clone/WA MRSA-60
CC5	CC5-MSSA
	CC5-MSSA [PVL+]
	ST228-MRSA-I, South German EMRSA/Italian Clone
	ST228-MRSA, South German EMRSA (subclone with truncated SCCmec element)
	ST5-MRSA-I, Geraldine Clone
	ST5/ST225-MRSA-II, Rhine-Hesse EMRSA
	ST5-MRSA-II [tst1/mer+], UK-EMRSA-3/Irish AR11
	ST5-MRSA-II [tst1+], UK-EMRSA-3/New York-Japan Clone
	ST5-MRSA-VII (SCC-JCSC6082)
	CC5/ST835-MRSA-(NovelSCCmec), WA MRSA-40/46
	CC5-MRSA-II [ACME+]
	CC5-MRSA-III
	CC5-MRSA-IV [ACME+]
	CC5-MRSA-IV&SCCfus
	CC5-MRSA-IV&SCCfus, "Maltese Clone"
CC5	CC5-MRSA-IV&VI
	CC5-MRSA-IV, Paediatric clone
	CC5-MRSA-IV, Paediatric clone [edinA+], WA MRSA-65
	CC5-MRSA-IV, Paediatric clone [PVL+/edinA+], WA MRSA-64
	CC5-MRSA-IV, Paediatric clone [PVL+]
	CC5-MRSA-IV, Paediatric clone [sed/j/r+]
	CC5-MRSA-IV, Paediatric clone [tst1+]
	CC5-MRSA-Ivar., WA MRSA-18/21/48/103
	CC5-MRSA-SCC(MRSAZH47)/IV&V
	CC5-MRSA-SCC(MRSAZH47)/IV&V [PVL+]
	CC5-MRSA-V [ACME+], WA MRSA-80
	CC5-MRSA-V [PVL+]
	CC5-MRSA-V [sec/d/j/l/r+], WA MRSA-87
	CC5-MRSA-V [sed/j/r+], WA MRSA-11/34/35/90
	CC5-MRSA-V [tst+]
	CC5-MRSA-V, WA MRSA-14/81/85/86
	CC5-MRSA-VI, New Paediatric Clone
	CC5-MRSA-VI, New Paediatric Clone [PVL+]
	CC5-MRSA with atypical SCCmec element "Montreal strain I"
	CC5-MRSA with atypical SCCmec element "Montreal strain II"



Clonal complex	Strain
CC6	CC6-MRSA-IV&V, WA MRSA-66
	CC6-MRSA-IV, WA MRSA-51
	CC6-MSSA
	CC6-MRSA-V
CC7	CC7-MSSA
	CC7-MRSA-IV
	CC7-MRSA-V
	CC7-MRSA-VI
CC7 (ST1048)	ST1048-MRSA-IV
CC8	CC8-MSSA
	CC8-MSSA [ACME+]
	CC8-MSSA [PVL+]
	CC8-MSSA-SCCfus
	ST247-MRSA-I, North German/Iberian EMRSA
	ST250-MRSA-I, Early/Ancstral MRSA
	ST8-MRSA-IIA/B/D, Irish AR13/14
	ST8-MRSA-IIA/B/D&SCC-M1, Irish AR13/14
	ST8-MRSA-IIC/E, Irish AR13/14
	ST8-MRSA-IIC/E&SCC-M1, Irish AR13/14
	CC8-MRSA-IV [sea+], Lyon Clone/UK-EMRSA-2
	CC8-MRSA-IV [sea+], UK-EMRSA-6
	CC8-MRSA-IV [tst1+]
	CC8-MRSA-IV&SCCfus
	CC8-MRSA-IV&V, WA MRSA-92
	CC8-MRSA-IV, [PVL+, sed/j/r+], WA MRSA-62
	CC8-MRSA-IV, Lyon Clone (sea-neg. Variant)/WA MRSA-88
	CC8-MRSA-IV, UK-EMRSA-14/WA MRSA-5
	CC8-MRSA-IV, USA500
	CC8-MRSA-V
	CC8-MRSA-V&SCCfus [ACME+], WA MRSA-77
	CC8-MRSA-V&VI
	CC8-MRSA-VIII
	ST254-MRSA-(atypical SCCmec), Hannover EMRSA (subclone)
	ST254-MRSA-IV&V, UK-EMRSA-10/Hannover EMRSA
	ST8-MRSA-(IVF+ccrA/B-4)/-VI
	ST8-MRSA-(IVG/E+ccrA/B-4), UK-EMRSA-12/13, Irish AR43
	ST8-MRSA-IV [PVL+/ACME-], ACME-negative Variant of USA300
	ST8-MRSA-IV [PVL+/ACME+], USA300
	ST8-MRSA-IV, putative PVL-deletion mutant of USA300
	CC8-MRSA-VI&SCCfus
CC8 (ST72)	ST72-MRSA-IV [PVL+], WA MRSA-44
	ST72-MRSA-IV, USA700
	ST72-MRSA-V, WA MRSA-91
	ST72-MSSA
CC8 (ST72)	ST72-MSSA [PVL+]
CC8 (ST239)	ST239-MSSA
	ST239-MRSA-III, Vienna/Hungarian/Brazilian Clone
	ST239-MRSA-III [ACME+]
CC9	CC9-MSSA
	CC9-MRSA-III
	CC9-MRSA-IV
	CC9-MRSA-V
	CC9-MRSA-V-atyp/truncated
	CC9-MRSA-IX
CC9 (ST834)	ST834-MSSA
	ST834-MRSA-IV [PVL+]
	ST834-MRSA-IV, WA MRSA-13
	ST834-MRSA-VI/SCCfus?
CC10	CC10-MSSA
CC12	CC12-MSSA
	CC12-MRSA-IV, WA MRSA-69
	CC12-MRSA, WA MRSA-59
CC15	CC15-MSSA [PVL+]
	CC15-MSSA

Clonal complex	Strain
	CC15-MRSA-I
CC20	CC20-MSSA
	CC20-MRSA-V
CC22	CC22-MSSA
	CC22-MSSA [PVL+]
	CC22-MSSA-SCCfus [PVL+]
	CC22-MRSA-V [PVL+]
	CC22-MRSA-IV, UK-EMRSA-15/Barnim EMRSA
	CC22-MRSA-IV [tst1+], UK-EMRSA-15/"Middle Eastern Variant"
	CC22-MRSA-IV [Q6GD50+], UK-EMRSA-15/"Maltese Variant"
	CC22-MRSA-IV [PVL+]
	CC22-MRSA-IV [ACME+], UK-EMRSA-15/"Dublin Variant"
	CC22-MRSA-IV [ACME+/PVL+]
	CC22-MRSA-V
CC25	CC25-MSSA
	CC25-MSSA [PVL+]
CC30	CC30-MSSA [PVL+]
	CC30-MSSA
	CC30-MRSA-IV&VI, WA MRSA-102
	ST36/39-MRSA-II, UK-EMRSA-16
	CC30-MRSA-IV [PVL+], Southwest Pacific Clone
	CC30-MRSA-IV [PVL-/tst1-]
	CC30-MRSA-IV [PVL-/tst1+], WA MRSA-68
	CC30-MRSA-IV&SCCfus
	CC30-MRSA-V
	CC30-MRSA-V [PVL+]
CC30 (ST34)	ST34-MSSA
CC45 (agr I)	CC45-MSSA
	CC45-MSSA [PVL+]
	CC45-MRSA-II, USA600
	CC45-MRSA-IV, Berlin EMRSA
	CC45-MRSA-IV [ACME+]
	CC45-MRSA-IV [tst1+/ACME+]
	CC45-MRSA-IV [tst1+]
	CC45-MRSA-IV&SCCfus
	CC45-MRSA-V
	CC45-MRSA-V [ACME+], WA MRSA-106
	CC45-MRSA-V [PVL+]
	CC45-MRSA-V [tst1+], WA MRSA-4
	CC45-MRSA-V&VI
CC45 (agr IV)	CC45/agrIV-MSSA
	CC45/agrIV-MRSA-IV, WA MRSA-23
	CC45/agrIV-MRSA-VT, WA MRSA-84
	CC45/agrIV-MRSA-IV&V
CC49	ST49-MSSA
	ST49-MRSA-V
	ST49-MSSA [PVL+]
CC50	CC50-MSSA
	CC50-MRSA-V&SCCfus
CC59	CC59-MSSA
	CC59-MSSA [PVL+]
	CC59-MRSA-IV [PVL+], USA1000
	ST59-MRSA-IV, WA MRSA-55/56 PVL-negative Variant
	ST59-MRSA-IV, WA MRSA-73
	ST59-MRSA-IV [PVL+], WA MRSA-55/56
	ST87-MRSA-IV, WA MRSA-24
	ST59/952-MRSA-V(T) [PVL+], Taiwan Clone
	ST59-MRSA-V
	CC59-MRSA-V [PVL+]
	ST59-MRSA-IV&V, WA MRSA-15
CC75 and related	CC75 lineage: MSSA, related to ST1223
	CC75 lineage: MSSA, related to ST1667
	CC75 lineage: MSSA, related to ST75
	ST1303-MRSA-IV, WA MRSA-76

Clonal complex	Strain
	ST1304-MRSA-IV, WA MRSA-72
	ST75-MRSA-IV, WA MRSA-8/79
	ST75-MRSA-V
	ST883-MRSA-IV, WA MRSA-47
CC80	CC80-MSSA [PVL+]
	CC80-MRSA-IV
	CC80-MRSA-IV [PVL+], European caMRSA Clone
	CC80-MRSA-(truncated/atypical SCCmec)
	CC80-MRSA-(truncated/atypical SCCmec) [PVL+]
CC88	CC88-MSSA
	CC88-MSSA [PVL+]
	CC88-MRSA-IV [etA+]
	CC88-MRSA-IV [PVL+]
	CC88-MRSA-IV, WA MRSA-2
	CC88-MRSA-V
	CC88-MRSA-V [PVL+]
	CC88-MRSA-VI
	CC88-MRSA-VII (SCC-JCSC6082)
	CC88-MRSA-Vtrunc. [PVL+]
ST93	ST93-MSSA
	ST93-MSSA [PVL+]
	ST93-MRSA-IV [PVL+], Queensland Clone
	ST93-MRSA-V
	ST93-MRSA-V [PVL+]
	ST93-MRSA-IV [PVL-]
CC96/154	CC96/154-MSSA
	CC96/154-MSSA [PVL+]
	ST154-MRSA-IV [PVL+], "Central Asian caMRSA"
CC97	CC97-MSSA
	CC97-MSSA-SCCmer
	CC97-MRSA-(I+V)
	CC97-MRSA-IV, WA MRSA-54/63
	CC97-MRSA-V
	CC97-MRSA-V&ACME
CC97 (ST71)	ST71-MSSA
CC101	CC101-MSSA
CC121	CC121-MSSA
	CC121-MSSA [PVL+]
	CC121-MRSA-IV
	CC121-MRSA-V [PVL+]
	CC121-MRSA-V, WA MRSA-22
	CC121-MRSA-VT, WA MRSA-93
CC130	CC130-MSSA/CC130-MRSA-XI
CC133	CC133-MSSA
CC133/ST2111	ST2111-MSSA
CC152	CC152-MSSA [PVL+]
	CC152-MRSA-V [PVL+]
CC182	CC182-MSSA [ccrA/B-2/kdp-]
	CC182-MSSA [ccrA/B-2/kdp+]
CC188	CC188-MRSA-IV, WA MRSA-38/78
	CC188-MRSA-V
	CC188-MSSA
	CC188-MSSA [PVL+]
CC188 (ST1774)	ST1774-MRSA-IV&ACME
CC361	CC361-MSSA
	CC361-MRSA-IV, WA MRSA-29
	CC361-MRSA-V, WA MRSA-70
	CC361-MRSA-VIII, WA MRSA-28
CC398	CC398-MSSA
	CC398-MSSA [PVL+]
	CC398-MRSA-IV
	CC398-MRSA-V
	CC398-MRSA-V [PVL+]
	CC398-MRSA-(truncated/atypical SCCmec)

Clonal complex	Strain
CC398 (ST291/813)	ST291/813-MSSA
	ST291/813-MSSA [PVL+]
CC479	CC479-MSSA
CC509	CC509/ST207-MRSA-V
CC509	CC509-MSSA
CC522	ST522-MSSA
CC692	CC692-MSSA
CC705 (=CC151)	CC705-MSSA
CC1021	CC1021-MSSA
	CC1021-MSSA [PVL+]
ST140	ST140-MRSA-IV
	ST140-MSSA
ST350	ST350-MSSA
ST425	ST425-MSSA/ST425-MRSA-XI
ST426	CC395-MSSA
	CC395-MRSA-IV
ST707	ST707-MSSA [ccrA/B-2/kdp-]
	ST707-MSSA [ccrA/B-2/kdp+]
ST779	ST779-MSSA
	ST779-MRSA-IV, WA MRSA-100
	ST779-MRSA-V
	ST779-MRSA-(novel SCCmec)
ST816	ST816-MSSA
ST913	CC913-MRSA-IV
ST942	ST942-MSSA
	ST942-MSSA [PVL+]
ST1093	ST1093-MSSA
ST1643	ST1643-MSSA
ST1755	ST1755-MSSA
ST1946	ST1946-MSSA/ST1946-MRSA-XI
Unknown /Singletons	ST(102-146-6-18-7-50-2)-MSSA
	ST(23-13-186-5-1-1-40)-MSSA
	ST(23-13-186-5-1-1-40)-MSSA [PVL+]